

Five Essential Tools to Achieve Surgical Convenience

HROSSEN
IMPLANT

Assist KIT

Five Essential Tools (KIT Components)


1. Tissue Punch 2. Bone Profiler 3. Trephine Drill 4. Crest Remover 5. Autobone Collector

1 Tissue Punch



Reduces Healing Time by Minimizing Gum Incisions

- Guide-pin configuration allows tissue removal at the exact location
- Tissue punch with a diameter 0.7-1.5 mm smaller than the healing abutment is recommended

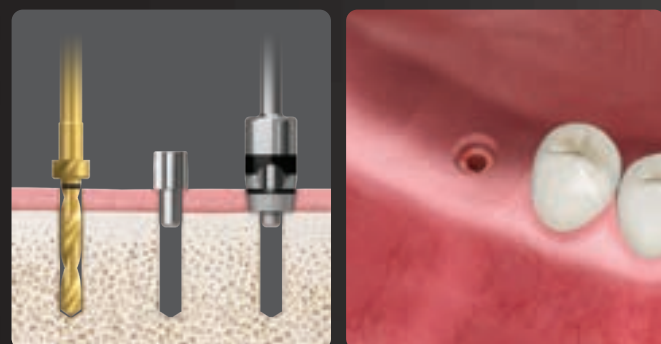


2 Bone Profiler




Convenient Bone Removal without Internal Damage of Implant

- Easier removal of protruding bones around the implant
- Detachable guide allows cleaning the bones without damaging the interior of the implant

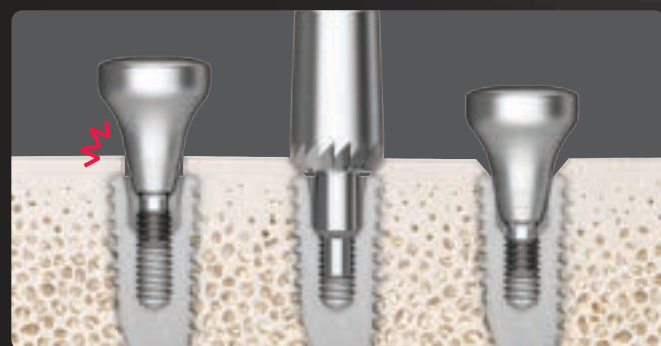


Bone Interference → Bone Removed → Normal Locking

3 Trephine Drill

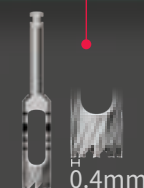

0.4 mm Blade Thickness to Minimize the Amount of Bone Removed

- Excellent cutting force allows easy acquirement of block bone or remove failed fixtures-0.4 mm thickness on one side minimizes the amount of bone removed (more than other companies)

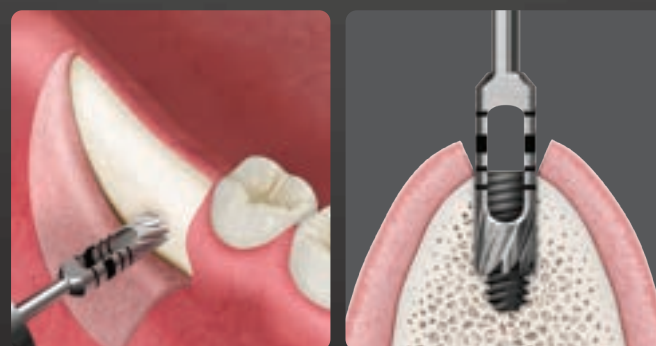


4 Crest Remover



Net-like Design Allows Fast and Simple Cleaning of the Ridge

- Blade is designed in a net shape to minimize the rim phenomenon compared to the round bur

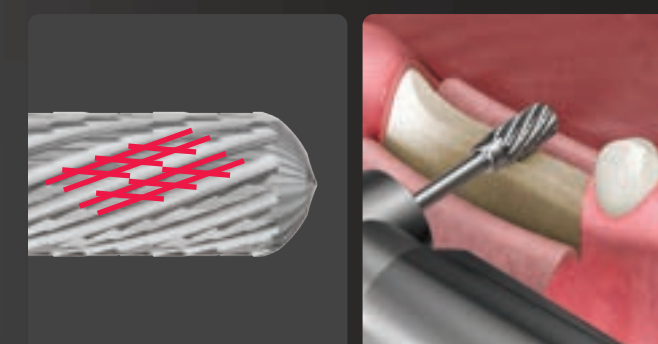

5 Autobone Collector



Easier Bone Collection Achieved Solely by Drilling

- Simple and quick to collect the chin bone with excellent cutting force and design
- Safer bone collecting procedure due to stopper (0.18 g acquired at one time)

